



The World's Most Portable  
Full-Function Dual-Channel  
Vascular Diagnostic System

# MULTILAB<sup>®</sup> SERIES 2-CP-R

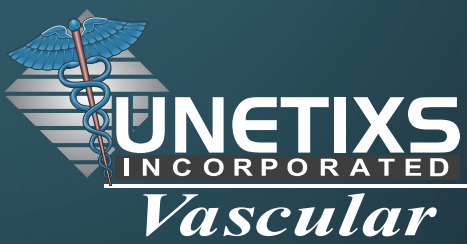

*Unetix, Inc.*

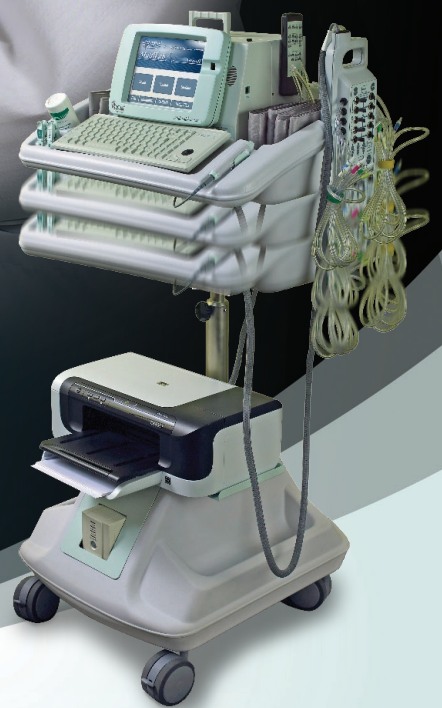



## Full-Function, Fully Portable

The MultiLab 2-CP-R combines capability with portability, offering the same technology as the acclaimed MultiLab LHS in an easy-to-carry unit. The 2-CP-R features TruDop® Doppler probes, calibrated PVR, dual-channel PPG probes with clips, and programmable software. You can generate a full array of lower and upper arterial studies, including ankle brachial index, toe brachial index, full segmental exams, post-exercise exams, full waveforms including Doppler, PVR and PPG, and venous reflux dual-channel exam.



## General Specifications

### Instrument Features:

- Composite Ergonomic TruDop® Doppler Probes
- Calibrated Dual-Channel Pulse-Volume Recording (C-PVR)
- Dual-Channel UltraSense Photoplethysmography (PPG)
- Soft-Touch Dual PPG Digit Clips
- 8.5" Diagonal T.F.T Color Touch Screen
- Full-Function Wireless Remote
- Weight-Complete: 18 lbs.
- Built-in CPU, Hard Disk Drive
- Three USB Ports
- Windows 10 Operating Platform
- Available DICOM
- Raynaud Study Capability
- Rear Panel Outputs for:
  - Color Printer, VGA Color Monitor
  - External Keyboard, Serial Port
  - Network-Compliant

### Optional Software:

- MultiLab® and Revo® Reporter Software
- Accreditation Software

Design and system specifications subject to change/enhancement without notice. System components and accessories may not appear as shown. Windows XP is a trademark of Microsoft Corporation.

### Standard Accessories:

- Color Page Printing (HP usually supplied)
- Multi-port Electronic Cuff Selector
- Full-Function External Remote Keyboard
- Complete Set of Pressure Cuffs to Perform All Studies
- Bilateral Testing Capabilities
- Custom Software-Configurable Protocols
- Custom Software-Configurable Test Reports

### Electronic Cuff Manifold Features:

- 10-Cuff Simultaneous Hookup
- Full-Function Remote Control for Automatic Inflation, Deflation, Capture, Forward, Back, Scrolling and more

### Cart Specifications:

- 21" x 24" x 40"
- 38 lbs. Total Weight
- 4" Swivel and Locking Casters
- Cuff and Probe Storage
- Miniature Full-Function Keyboard
- Printer Compartment
- Ultrasound Gel Holder
- Integral Cuff Manifold Holder

(800) 486-3849  
www.unetixs.com

Unetix Vascular, Inc.  
125 Commerce Park Road  
North Kingstown, RI 02852  
401-294-7559 Tel  
401-294-3893 Fax



Comprehensive step-by-step reference manual for performing noninvasive peripheral vascular procedures. Includes 5.25 CMEs.

**U.S. Inquiries:** sales@unetixs.com  
**International Inquiries:** international@unetixs.com

**CSJ CRITICARE**  
SYSTEMS, INC.

SM 165 Rev A